

Individually matched to your specific requirements.

The special feature of the G 2.5 TTL Galilean system is that both the eyepiece lens and the carrier lens are available in prescription form. Your expert ZEISS optician will be happy to advise you on the combination which is perfectly adapted to your specific working requirements.

Relaxed work without strain.

ZEISS prism loupes and head-mounted systems enable you to work in comfort and offer outstanding perception of detail. The ideal relationship between magnification and working distance ensures a relaxed, almost upright posture – for every body size in any working situation. The result is an individual solution which offers you ideal support in your daily work.

Crafted to your individual centring data.

All G 2.5 TTL teloupe systems are unique and feature the highest ZEISS quality standards. Your ZEISS optician measures your individual centring data and, like all ZEISS special visual devices, the system is then produced especially for you. Benefit now from this vision aids - a wealth of individuality which offers many benefits.


EN_20_010_60711

Carl Zeiss Vision GmbH
www.vision.zeiss.com
info-de@vision.zeiss.com

See more. Live more.
ZEISS precision loupes.



A leader in precision optics since 1846.

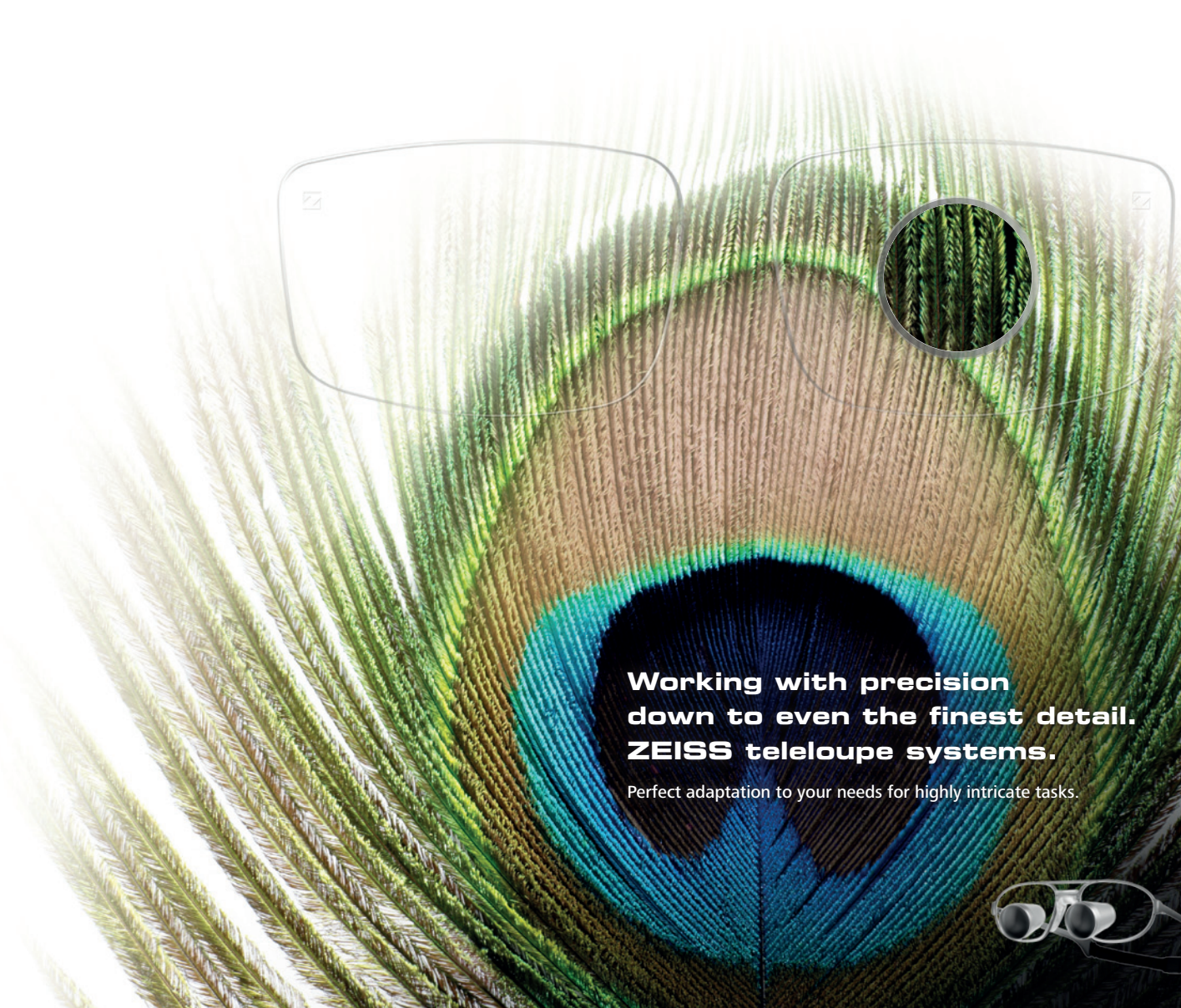


**Working with precision
down to even the finest detail.
ZEISS teloupe systems.**

Perfect adaptation to your needs for highly intricate tasks.

**Every detail matters.
With ZEISS precision loupes
you won't miss a thing.**

Great news for all those whose work demands an eye for even the finest detail: the new ZEISS G 2.5 TTL teleloupe system is now available with all the technical features required for optimum visual performance – whatever the task and whatever the situation.

Outstanding image quality over a large field of view – the typical feature of all magnifying visual devices from ZEISS. The latest milestone: the ZEISS G 2.5 TTL teleloupe system with a field of view which is 30% larger than the GTX system. **Make your everyday working life easier with the ZEISS G 2.5 TTL teleloupe system.**

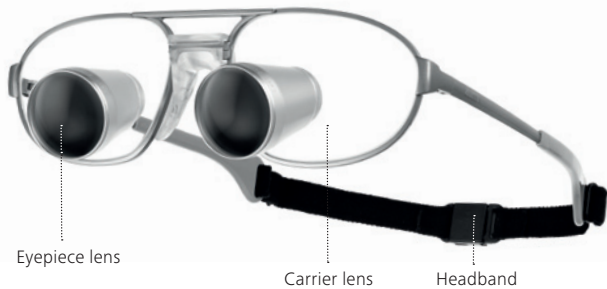

Technical data:

- Magnification: 2.5x
- Working distances: 350 / 400 / 450 mm
- Fields of view: 75 / 85 / 95 mm



Precision in every detail

ZEISS G 2.5 TTL with 2.5x magnification and Saphiro² LED illumination



G 2.5 TTL: the benefits in your daily work

Comfort

- System mount with straight sides and headband enables perfect adjustment to the head
- Fatigue-free work even over longer periods

Performance

- Optimum overview of the work area
- Distortion-free and high-contrast image
- Ergonomic working conditions thanks to ideal matching of magnification and working distances

Individuality

- All individual centring data taken into consideration
- Made-to-measure fitting of the system into the carrier lenses
- Easy to use without any need to make adjustments
- Available in the following product versions:
 - TTL basic:** the eyepieces and the carrier spectacles are fitted with normal lenses
 - TTL pro 1:** the eyepieces are fitted with prescription lenses
 - TTL pro 2:** both the eyepieces and the carrier spectacles are fitted with prescription lenses
- Coaxial illumination as an option